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NOTIZEN

DEAR CUSTOMER,

Congratulations on your new hearing systems. You have chosen sophisticated, high-quality products. By purchasing these hearing systems you have taken an important step to enable yourself to hear well again. We are sure these hearing systems will soon become indispensable to you and make your life safer and more enjoyable. This user manual will help you to familiarise yourself with your new hearing system.

If you require further information, please contact your hearing system specialist who will be happy to help.

THE HEARING SYSTEMS	36-37
INTENDED USE	38
SAFETY INSTRUCTIONS	39-41
HANDLING THE BATTERIES	42-44
– Battery type	42
– Fading battery / rechargeable	43
– Inserting the battery / rechargeable	43
– Removing the battery / rechargeable	44
INSERTING AND REMOVING THE HEARING SYSTEM	45
– Inserting the hearing system	45
– Removing the hearing system	45
SWITCHING ON AND OFF	46-47
– With the battery compartment	46
– With program button (Mood)	46
– With sliding switch (Sun)	46
– With the remote control	47
– Power-on delay	47
PROGRAM SELECTION	48
– With program button (Mood)	48
– With sliding switch (Sun)	48
– With the remote control	48
TELEPHONE INDUCTOR AND MIXED PROGRAM (SUN)	49

CONTROLLING THE VOLUME	50
– Automatic	50
– With sliding switch (Sun)	50
– With the remote control	50
WIRELESS	51
– Wireless functionality	51
– Technical data for wireless hearing systems	51
MAINTENANCE AND CARE	52-54
– Standard ear moulds (domes) and customised ear moulds (CLIC MOULDS)	53
– Cleaning ear moulds	54
ACCESSORIES	55-56
– Battery charger station (for selected hearing system models)	55
– Bluetooth remote control (for selected hearing system models)	55
– Remote control (optional)	56
TROUBLESHOOTING	57
DISPOSAL / CE	58
– Disposal	58
– CE-mark	58
CONFIGURATION OF THE HEARING SYSTEMS	59-60

THE HEARING SYSTEMS

MOOD

Battery type 312
rechargeable

RIXX

Battery type 10
not rechargeable

SUN

Battery type 13
rechargeable

Please tick customised hearing system.

1	Main hearing system component
2	Receiver
3	Microphone protector
4	Battery compartment
5	Charge contact (not for all models)
6	CLIC MOULD (customised earpiece as alternative to the dome)
6a	Removal line
7	Receiver connector
8	Dome (standard earpiece)
9	Program button
10	Sliding switch Can be programmed by your hearing system specialist as: <ul style="list-style-type: none"> – Program button – Program button with on/off function (standard) – Volume control

INTENDED USE

Hearing systems are intended to improve the hearing of hearing-impaired persons. The diagnosis and prescription of hearing systems must be performed by a qualified person, e.g. by ENT specialists, audiologists or hearing system specialists.

Please read and follow the instructions in this user manual to prevent injuries and protect your hearing systems from damage.

Symbols:

Take note of the information accompanied by the **WARNING**, **CAUTION** or **NOTE** warning symbols.

WARNING draws attention to a situation in which there is a risk of serious injury.
CAUTION draws attention to a situation in which there is a risk of medium or minor injury.

NOTE

NOTE informs the user of a possible risk of damage.

Recommendations and tips for an improved experience with your hearing system.

Instruction: indicates that an action must be performed.

Indicates an option that is not available on every hearing system. Ask your hearing system specialist to mark the relevant information.

SAFETY INSTRUCTIONS

NOTE

Your hearing system is sensitive to extreme heat, high air humidity, strong magnetic fields (> 0.1 T), X-rays and mechanical stress.

- Do not expose your hearing system to extreme temperatures or significant moisture.
- Do not leave it lying in the sun.
- Do not wear it when in the shower, or when putting on your make-up, perfume or aftershave, hairspray or sun-tan lotion.
- Do not wear your hearing systems when undergoing microwave therapy or medical procedures in which strong magnetic fields, RF fields or X-rays occur.
- Do not put your hearing systems in the microwave.

**WARNING**

Risk of choking if you swallow small parts.

- Keep your hearing systems, batteries and accessories away from children.
- If parts are swallowed, consult a doctor or go to the hospital immediately.
- This hearing instrument is not designed for children under the age of 36 months or mentally disabled persons. For suitable products please refer to your Hearing Care Professional.

**WARNING**

Danger of explosion!

- Do not wear your hearing systems in environments in which there is a risk of explosion (e. g. in a mine).

NOTE

Leaking batteries / rechargeables will damage the hearing systems.

- Switch the hearing systems off when you are not wearing them to preserve the battery / the rechargeable.
- Remove the batteries / rechargeables if you will not be wearing the hearing systems for some time.

NOTE

Batteries contain hazardous, polluting substances.

- Do not dispose of old batteries in the domestic waste.
- Dispose of the batteries in accordance with current regulations, or hand them to your hearing system specialist.

- In some countries, the use of radio devices is subject to certain restrictions.
- Further information may be obtained from your local authority.

HANDLING THE BATTERIES / RECHARGEABLES

Original size of rechargeable/
not rechargeable battery type

i Use only batteries / rechargeables of the correct battery size.

Charge a rechargeable battery before using it for the first time. (Further information about this may be found in the battery charger operating instructions.)

Do not remove the protective film from the batteries until you are ready to use them. The zinc-air battery then needs a few minutes to reach full power.

The positive side is always marked. You will recognise the negative side by the ring let into it.

Always keep spare batteries / rechargeables to hand.

FADING BATTERY / RECHARGEABLE


Your hearing system specialist can program your hearing system so that you are warned by an acoustic signal as the battery charge / the rechargeable charge declines. Change the battery or charge the rechargeable when you hear this acoustic signal or if the amplification is decreasing.

INSERTING THE BATTERY / THE RECHARGEABLE

Open the battery compartment fully with a finger nail (see illustration).

Insert the correct battery / rechargeable. The plus sign on the battery faces up (see illustration).

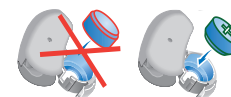Close the battery compartment. If you cannot close the battery compartment easily, check that the battery / rechargeable has been inserted properly. Please do not use any force when closing the battery compartment. The hearing system is ready to use when the battery compartment is fully closed.

REMOVING THE BATTERY / THE RECHARGEABLE


- Open the battery compartment fully.
- Take care not to damage the battery compartment as the battery is removed.

INSERTING AND REMOVING THE HEARING SYSTEM


INSERTING THE HEARING SYSTEM

- Switch your hearing system off.
- Grip the ear mould between thumb and index finger and insert it carefully into the auditory canal. It may help if you pull your ear lightly up or backwards.
- Fit the hearing system behind the ear.
- Switch the hearing system on.

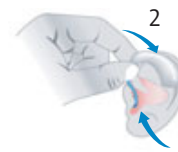REMOVING THE HEARING SYSTEM

- Grip the ear mould between thumb and index finger and pull it carefully out of the ear. It may help if you push your ear gently from rear to the front.

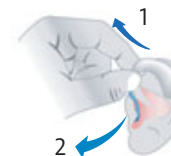i Your hearing system specialist will affix colored markings in the battery compartment to identify the sides. A red marking in your battery compartment identifies the right-hand hearing system, with a blue marking identifying the left-hand hearing system.

SWITCHING ON AND OFF

Ask your hearing system specialist to mark the options relevant for you.

WITH THE BATTERY COMPARTMENT

ON → Close the battery compartment.

OFF → Open the battery compartment.

WITH PROGRAM BUTTON (MOOD)

ON → Press the program button and hold for 2 seconds.

OFF → Press the program button and hold for 2 seconds.

WITH SLIDING SWITCH (SUN)

ON → Slide the switch upwards or downwards for 2 seconds.

OFF → Slide the switch upwards or downwards for 2 seconds.

WITH THE REMOTE CONTROL

See the remote control user manual.

POWER-ON DELAY

Your hearing system has a power-on delay which your hearing system specialist is able to activate. This ensures that the hearing system only switches itself on after an automatic delay of approximately 12 seconds to prevent feedback whistling as you insert the hearing system in the ear.

PROGRAM SELECTION

WITH PROGRAM BUTTON (MOOD)

Press the program button. The hearing system changes to the next program.

WITH SLIDING SWITCH (SUN)

Slide the switch upwards once to change to the next program.

Slide the switch downwards once to change to the previous program.

If the sliding switch is configured to fulfill the “program selection” function, it is not possible to use it for controlling the volume.

WITH THE REMOTE CONTROL

See the remote control user manual.

TELEPHONE INDUCTOR AND MIXED PROGRAM (SUN)

Your hearing system Sun is equipped with a telephone inductor, which makes it easier for you to telephone. The telephone inductor can also be used for better understanding in buildings in which induction loops are installed (e.g. church, theatre, concert hall, cinema, etc.).

If you want to telephone and at the same time understand what someone nearby is saying, you can change to a mixed program, in which both microphone and telephone inductor are activated.

ATTENTION

Not all telephones are suitable for use with a hearing system. This applies in particular to cordless telephones and cell phones.

Using the sliding switch on top of the hearing system you can select the hearing program that is suitable for the particular situation, or change from the microphone program to the telephone or mixed program.

The change of program is accompanied by an audible signal consisting of up to 5 successive beeps (depending on the model).

Please ask your hearing system specialist about this.

CONTROLLING THE **VOLUME**

The volume of your hearing systems is controlled automatically. It is normally not necessary for you to modify the volume manually. If you prefer manual operation, your hearing system specialist is able to change your hearing systems over to manual operation.

AUTOMATIC

The volume is set automatically when the hearing systems are used.

USING SLIDING SWITCH (SUN)

Slide the switch upwards to increase the volume.

Slide the switch downwards to decrease the volume.

If the sliding switch is configured to fulfill the 'volume control' function, it is not possible to use it for changing hearing programs.

WITH THE REMOTE CONTROL

See the remote control user manual.

WIRELESS FUNCTIONALITY

Your hearing system has wireless functionality.

- Permits the use of a remote control.
- Automatically adjusts both hearing systems at the same time if you change the system settings, e. g. volume or hearing program.

i The effectiveness of the wireless connection can be impaired by electromagnetic radiation, e. g. from computer screens.

- Increase your distance from the source of interference if you suffer interference.

TECHNICAL DATA FOR WIRELESS HEARING SYSTEMS

Frequency range	3.28 MHz
Pulse duration	325 ms
Transmitting power	42 µW

MAINTENANCE AND CARE

Your hearing system is highly sensitive. Optimum performance can be ensured by proper handling and care. Clean your hearing system regularly to prevent damage to the system and health problems.

NOTE

Never hold your hearing systems under water.

- Use a soft, dry cloth to clean the hearing system, as necessary.
- If your hearing system is distorting or sounds weak, have your hearing system specialist replace the microphone cover (upper housing of the hearing system).
- Store your hearing system in a dry bag overnight. This protects the systems against harmful humidity.
- Ask your hearing system specialist for special care kits and find out what you can do to keep the hearing system in good condition.

STANDARD EAR MOULDS (DOMES) AND CUSTOMISED EAR MOULDS (CLIC MOULDS)

CAUTION

Risk of injury from using external receiver and the ear mould incorrectly.

- Always wear the external receiver with an ear mould.
- Make sure that the ear mould is properly connected to the external receiver.
- If your ear mould becomes stuck in your ear, have it removed by a medical expert.
- Fit the external receiver and the ear moulds carefully and not too deeply into the ear.
- Wear the external receiver and the ear mould only as described in this manual.

Your hearing system has one of the following ear moulds:

Click Dome

CLIC MOULD


CLEANING EAR MOULDS

Cerumen (ear wax) on the ear mould degrades the sound quality.

Always clean the ear moulds immediately after you have removed the hearing systems from your ear. The cerumen will then not dry out and harden.

Remove the cerumen with a soft cloth. No other tools are required.

If the cerumen has dried on, moisten a soft cloth with cerumen cleaning fluid. Dab the ear mould with the cloth and wait a moment until the cerumen is loosened.

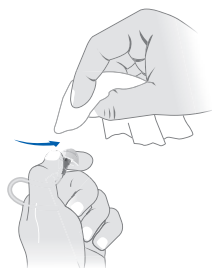ACCESSORIES

BATTERY CHARGER STATION (FOR SELECTED HEARING SYSTEM MODELS)

Your hearing systems work with normal batteries or rechargeables. You can charge the batteries for one or two hearing systems in this easy-to-use battery charging station. The charging station switches itself off automatically after charging for six hours and the hearing systems are then ready for use for another day.

You can get further information from your hearing system specialist or in the user manual for the battery charging station.

BLUETOOTH REMOTE CONTROL (FOR SELECTED HEARING SYSTEM MODELS)

A special Bluetooth connection makes it possible for selected hearing systems to receive audio signals from your mobile phone, television or MP3 player.

Further information is available from your hearing system specialist.

REMOTE CONTROL (OPTIONAL)

The wireless function in your hearing systems makes it possible for you to change the program and adjust the volume using the remote control.

You can get further information from your hearing system specialist or in the instruction manual for the remote control.

i Please consult your hearing system specialist if your hearing systems are being interfered with by someone else's remote control.

TROUBLESHOOTING

Problem	Solution
The sound is weak.	Replace the flat battery or charge the rechargeable. Clean the ear mould. Have the microphone cover replaced. Increase the volume.
The hearing system is whistling.	Insert the ear mould so that it is firmly and securely seated. Reduce the volume. Clean the ear mould.
Sound is distorted.	Reduce the volume. Replace the flat battery or recharge the rechargeable. Clean the ear mould.
You can hear signals from the hearing system.	Replace the flat battery or recharge the rechargeable.
Hearing system is not working.	Switch the hearing system on. Carefully close the battery compartment fully. Replace the flat battery or recharge the rechargeable. Make sure that the battery / rechargeable has been inserted correctly. Power-on delay is active. Wait approximately 12 seconds and then check again.

Consult your hearing system specialist if you have any further questions.

DISPOSAL

Within the European Union, the accessories listed are covered by Directive 2002/96/EC of the European Parliament and of the Council dated 27 January 2003 on waste electrical and electronic equipment. Most recently modified by Directive 2003/108/EC (WEEE).

Dispose of hearing systems, batteries and accessories in accordance with the regulations.

CE-MARK

With the CE-mark, Audio Service confirms compliance with European Directive 93/42/EEC on medical devices.

For products with Wireless 2.0, Audio Service confirms compliance with European Directive 99/5/EC (R&TTE) concerning radio equipment and telecommunications terminal equipment.

Subject to technical changes and errors.

CONFIGURATION OF THE HEARING SYSTEMS

Hearing system			
left			right
		Model	
		Serial number	
		Battery size	
<input type="checkbox"/> yes <input type="checkbox"/> no		Power-on delay	<input type="checkbox"/> yes <input type="checkbox"/> no
<input type="checkbox"/> yes <input type="checkbox"/> no		Program button	<input type="checkbox"/> yes <input type="checkbox"/> no
		Sliding switch as:	
<input type="checkbox"/> yes <input type="checkbox"/> no		Program button	<input type="checkbox"/> yes <input type="checkbox"/> no
<input type="checkbox"/> yes <input type="checkbox"/> no		Program button with on/off function	<input type="checkbox"/> yes <input type="checkbox"/> no
<input type="checkbox"/> yes <input type="checkbox"/> no		Volume control	<input type="checkbox"/> yes <input type="checkbox"/> no
Signal tones:			
<input type="checkbox"/> on <input type="checkbox"/> off		Program change	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off		Battery warning signal	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off		Volume limit	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off		Volume change	<input type="checkbox"/> on <input type="checkbox"/> off

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Audio Service GmbH · Zeppelinstraße 9 · D-32051 Herford